


PERCENT POTASSIUM CHLORATE IN A MIXTURE

Once we are able to describe pure substances in terms of atoms and molecules, and in terms of compounds and elements, we can begin to understand chemical reactions in a quantitative manner. We can predict product amounts and develop economical ways to produce the many synthetic substances we need. Furthermore, we have the ability to analyze substances based upon the amount of products formed in a reaction. This experiment is an illustration of just such a technique.

Potassium chlorate (KClO_3) decomposes on heating to produce potassium chloride and oxygen according to the equation below:

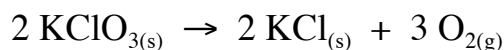

Since the decomposition of 2 mole of KClO_3 releases 3 moles of oxygen, you can determine the amount of KClO_3 in a sample from the weight loss due to loss of oxygen. In this experiment, you will determine the percent KClO_3 in an unknown mixture containing an unreactive chloride. By converting the weight loss to moles of oxygen lost, we calculate the moles of potassium chlorate in the starting material. From the moles of potassium chlorate we can convert to mass and then the percent of potassium chlorate in the unknown mixture.

A catalyst, MnO_2 (manganese (IV) oxide), is added to the reaction mix to speed up the reaction. Like all catalysts, the same amount of catalyst is present at the end of the reaction as in the beginning. Therefore, we will include the mass of the catalyst in with the mass of the test tube.

In a thermal decomposition such as this, the product must be heated to a constant weight before you can be sure the decomposition is complete. After the first heating, cooling and weighing, the sample must be heated again, cooled and reweighed. This process is continued until two successive weighing are within 5 mg of each other (three or four heatings.)

PROCEDURE: - *You must wear safety goggles while performing this lab!*

1. Set up a ring stand with a triangle as demonstrated by your teacher. The small white crucible should fit inside the triangle.
2. Place about 0.5 g of manganese(IV) oxide into a clean, dry small white crucible. Heat the crucible and catalyst with a Bunsen burner for about 3 minutes to drive off any moisture that may be in the catalyst or crucible.
3. When the crucible is cool enough to touch, record the entire mass to the nearest 0.001g. (The mass of a hot object will be greater than the mass of a cool object so please make sure that your crucible is cool before determining the mass.)
4. Add between 2.0 to 2.5 grams of the unknown mixture to the crucible. Mix the contents to obtain a somewhat uniform mixture. Record the mass of the crucible plus catalyst plus mixture to the nearest 0.001 g. Be sure to also record your unknown number!
5. Begin heating the crucible gently at first followed by a more aggressive treatment for a total of 10 minutes. Be aware that the sample may begin to bubble and spurt; if this happens, turn the flame down a bit.
6. Allow the sample to cool to room temperature. Record the mass to the nearest 0.001 g.
7. Reheat your sample for 5 minutes. Cool and record the mass. If your mass is within 0.005 g of the mass after the first heating, congratulations, you can move on to calculations. If not, you must reheat, cool, and weigh until you have two successive masses within 0.005 g of each other.
8. Clean up! Complete postlab questions.

PERCENT POTASSIUM CHLORATE

Name: _____

Lab Partner: _____

DATA:

Unknown Number _____

1. Mass of crucible + catalyst
(after drying moisture & before adding unknown) _____
2. Mass of crucible + catalyst + unknown
(before heating) _____
3. Mass after first heating _____
4. Mass after second heating _____
Mass after third heating (if necessary) _____
Mass after fourth heating (if necessary) _____

CALCULATIONS: (Show set up for each calculation)

7. Mass of oxygen lost _____
8. Mass of potassium
chlorate decomposed _____
9. Mass of unknown _____
10. % KClO_3 in unknown mixture _____

